

Patent-Pending

SYMKnife®

Multi
Grid
Holder


**What's
this?**

Diamond knife by SYNTEK, Japan

Benefits of using "MULTI-GRID HOLDER"


- 1) It will make your works much easier than now
- 2) Your hand will be released finally
- 3) For everybody from beginner to experts

How to use Multi-Grid Holder

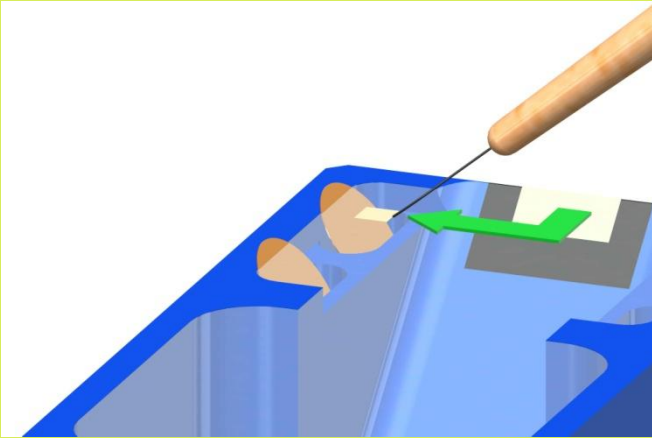


Step1

Put the grids into holders (up to 4).
Cut sample is guided to each grid.

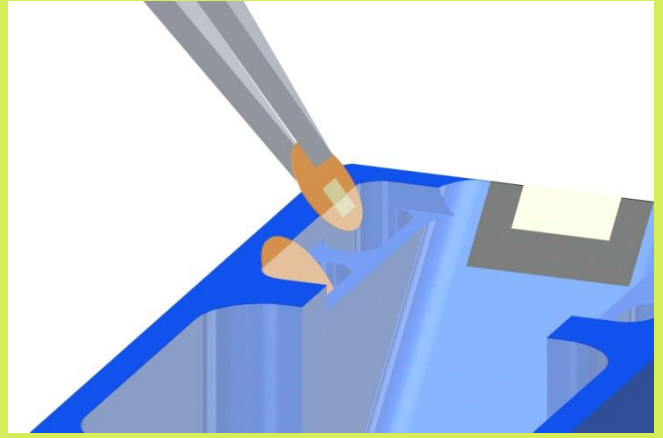


Step2

Collect first grid.

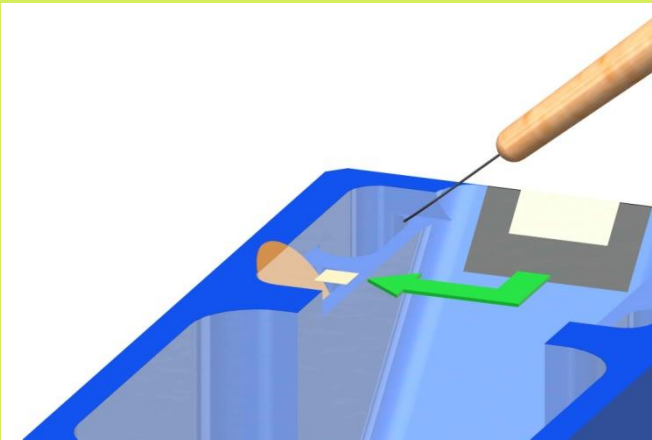

Step3

Cut sample is guided to another grid.

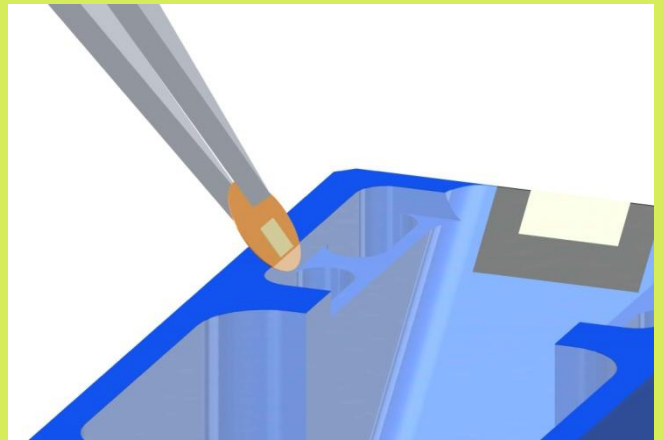

Step4

Collect second grid.

Description	Angle	Width(mm)	Model No.	New	Exchange
Ultra (Wet) The recommended range 200nm to 30nm	35°	2.0 (2.0-2.2)	MGH2035	○	○
	45°	2.0 (2.0-2.2)	MGH2045	○	○
		2.5 (2.5-2.7)	MGH2545	○	○
		3.0 (3.0-3.2)	MGH3045	○	○

*Replacement is only the same specification (diamond width) as yours.

*For improvement, some changes in appearance or specifications may be made.

*For the preparation of 30-40nm specimens, it varies depending on conditions such as sample width, embedding condition, hardness and so on.



SYNTEK CO., LTD.

info@syntek.co.jp

<http://www.syntek.co.jp/en>